

## B/L SACC (Supra Genual Cingulate Cortex )


Voxel Size RL X AP X SI 20 X 30 X 15


- Reformate Axial and Coronal Images from Sag T1 BRAVO.
- Axial and Coronal images are parallel to AC-PC line and RL/AP/SI tilt corrected.
- Make an Oblique axial slice parallel to Supra Genual Cingulate Cortex (SACC)








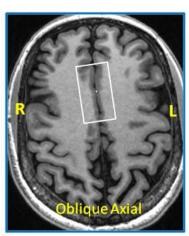

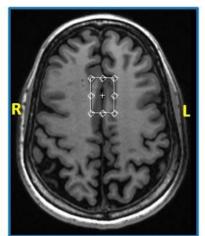



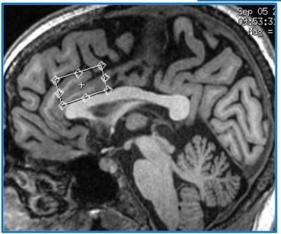

## Exam Rx:


- -On oblique axial place the voxel centered on the dot/cursor.
- Open all the three planes Axial, Coronal and Sag.
- Make sure B/L SACC is covered on all three planes
- -Stay away from corpus callosum, include grey matter of SACC.






## LT.DLPFC 25 X 25 X 15 (AP/SI/RL)


✓ Correct Tilts: R/L - S/I - A/P✓ Oblique to AC-PC

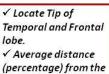



to the posterior vertical boundary of the DLPFC is 20%. ✓ DLPFC is located between Superior and Inferior frontal Sulcus. √ Place the voxel on the double oblique image √ Stay away from


Bone

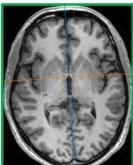


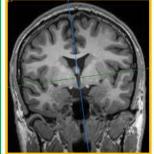


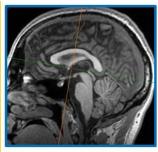


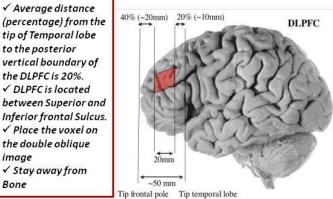


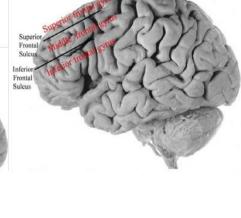


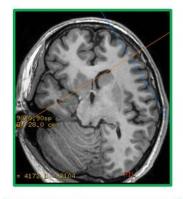

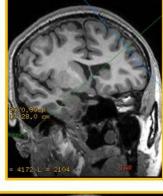

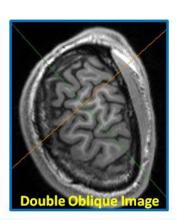

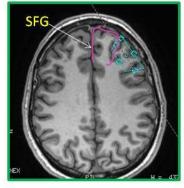

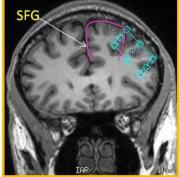



Voxel have to be placed inferior to Superior Frontal Gyrus (SFG)